

BENCH PROTECTORS

Clean Sheets Bench, Shelf, and Drawer Liner. White cross-linked polyethylene foam, 1/8" (3 mm) thick material. Comes in 50 foot (15 m) rolls in two widths; fits drawers (18 in) and shelves (12 in). Use it to line sinks, drawers, benches, shelves – anywhere protection is needed against chemicals, vibration, heat, dirt, and breakage. Can be cut to any size and stapled, tacked or taped.



Width	NCL #	Price/ea
12"	BN-628A	134.50
18"	BN-628B	186.50



DISPOSABLE BENCH LINER

Disposable Labmat Bench and Drawer Liner simplifies clean-up chores and safeguards work surfaces from stains, spills, and wear. The flexible polyethylene film backing is non-skid and prevents soak-through of most solutions. The highly absorbent upper surface traps water, dyes, corrosive chemicals, and radioactive solutions. The entire sheet can be rolled up and discarded after contamination. Labmat Liner is 20" wide and is packed in 50 ft rolls, each in a handy dispenser carton. Required footage is easily cut off with scissors while the roll is retained in the dispenser.


BL-101

50.00/rl

LABORATORY SAFETY MAT

This grooved laboratory safety mat is made of non-slip NONSLIP silicon rubber. It is temperature resistant up to 260°C.

Size	Color	NCL #	Price/ea
10" x 10"	Red	MR-250	15.95
10" x 10"	Black	MB-250	15.95
14" x 14"	Red	MR-350	21.60
14" x 14"	Black	MB-350	21.60

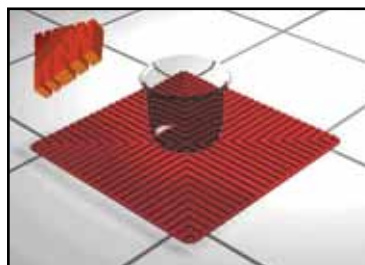



ACRYLIC LEVELING TABLE

Keep liquids perfectly level with this heavy-duty leveling surface. Each table has a leveling bubble and three leveling feet.

Size	NCL #	Price/ea
8" x 8"	LT-88	79.00
8" x 12"	LT-812	89.00

CERAMIC BENCH PROTECTORS

Protect expensive countertops from high temperatures associated with crucibles, glassware, etc. Each "thermoplate" is manufactured from glass ceramic and has corner feet to isolate it from the countertop. Stable over a wide temperature range (-392° to 1200°F), near zero thermal expansion and insensitive to thermal shock. Dark brown in color, these are stain-resistant and easy to clean. Well worth the money.



Size	NCL #	Price/ea
9" x 9"	ET-103	57.55
12" x 12"	ET-104	74.65



DRAWER LINER

18" x 18" soft, ribbed material fits in drawers or on shelves and keeps fragile items from rolling when drawers are opened.

DL-144 20.00/ea