

WTW FIELD DO/pH METER

Portable waterproof pH/DO/mv/Temp meter developed for hard on-site work. Features a 120-point datalogger, large LCD display, and RS232 interface. pH mode features automatic temperature compensation, and autorecognition of 4, 7, and 10 buffers, using a 2-point pH calibration. DO probe can be easily air-calibrated. Battery operated using four disposable AA batteries. Datalogging intervals are between 5 seconds and one hour. Logged data can be viewed from meter display or can be downloaded to a PC using optional software and cable. Both probes have 10-foot cables and are submersible. 3-Year warranty on meter.

**Includes:**

- ◆ Meter
- ◆ Sentex 41 pH probe with 10 ft cable
- ◆ CelloX[®] 325 DO probe with 10 ft cable
- ◆ pH 4, 7, and 10 buffers
- ◆ pH storage solution
- ◆ 3 membranes caps for DO probe
- ◆ Electrolyte filling solution for DO probe
- ◆ Cleaning solution for DO probe
- ◆ Calibration and storage sleeve for DO probe
- ◆ Polishing paper for DO probe
- ◆ Carrying case



Modes	Ranges	
DO (mg/l) & Temp	0.00 to 19.99	0 to 40°C
DO (% sat.) & Temp	0.00% to 99.99%	0 to 40°C
pH & Temp	-2.0 to +16.00	-5°C to 99.9°C
mv & Temp	-1250 to +1250	-5°C to 99.9°C

Description	NCL #	Price
Complete kit, as described above	WT-410	2090.00
Replacement Sentex 41 pH/Temp probe with 10 ft cable	WT-636	272.00
CelloX [®] 325 DO probe with 10 ft cable	WT-545	787.00
Replacement membrane caps for DO probe (3/pk)	WT-725	115.00
Electrolyte filling solution for DO probe (50 ml bottle)	WT-217	55.00
Cleaning solution for DO probe (50 ml bottle)	WT-204	52.00
Polishing paper for DO probe	WT-680	25.00
DO service kit, includes: 3 membrane caps, filling solution, cleaning solution, and polishing paper	WT-706	188.00



ORION® 4-STAR PORTABLE DO/pH METER

- ★ Measures DO, pH, and Temp
- ★ Bright LCD display with ready indication
- ★ IP67 waterproof
- ★ Autoread/hold function
- ★ Advanced self-tests and self-diagnostics on meter
- ★ 200 point datalogger which is downloadable to printer or computer via an optional RS232 or USB interface accessory
- ★ Three-year replacement warranty



DO:

- ★ Automatic barometric pressure compensation
- ★ Can be calibrated to Winkler value
- ★ Stability indicator
- ★ Salinity correction

pH:

- ★ Automatic or manual calibration (1 to 5 points)
- ★ Automatic buffer recognition

ORION® 4-STAR PORTABLE pH/DO METER ONLY

DP-1216 868.50/ea

PORTABLE pH/DO METER KIT II

Same as above, except different probes (with 4 ft cables): 083005MD DO probe and 9017BNMD pH/ATC triode.

DP-1216A 1339.20/ea

Specifications

pH	Range	-2.000 to 19.999
	Resolution	0.1 / 0.01 / 0.001
	Relative Accuracy	±0.002
DO	Concentration	0.00 to 90.00 mg/l
	Relative Accuracy	0.2 mg/l
	% Saturation	0.0 to 600
	Salinity Factor	0 to 45 ppt
	Barometric Pressure	450-850 mm Hg
	Temperature	-5 to 105°C

ORION® 4-STAR BENCHTOP DO/pH METER



- ★ Measures DO, pH, and Temp
- ★ Bright LCD display with ready indication
- ★ IP54 dustproof and splashproof
- ★ Autoread/hold function
- ★ Advanced self-tests and self-diagnostics on meter
- ★ 200 point datalogger which is downloadable to printer or computer via an optional RS232 or USB interface accessory
- ★ Three-year replacement warranty.



DO:

- ★ Automatic barometric pressure compensation
- ★ Can be calibrated to Winkler value
- ★ Stability indicator
- ★ Salinity correction

pH:

- ★ Automatic or manual calibration (1 to 5 points)
- ★ Automatic buffer recognition

Specifications

pH	Range	-2.000 to 19.999
	Resolution	0.1 / 0.01 / 0.001
	Relative Accuracy	±0.002
DO	Concentration	0.00 to 90.00 mg/l
	Relative Accuracy	0.2 mg/l
	% Saturation	0.0 to 600
	Salinity Factor	0 to 45 ppt
	Barometric Pressure	450-850 mm Hg
	Temperature	-5 to 105°C

ORION® 4-STAR pH/DO BENCHTOP METER ONLY

DP-1116 1089.00/ea

BENCHTOP pH/DO KIT

Includes: meter, 8102BNU ROSS Ultra® pH Electrode, 081010MD epoxy DO probe with BOD funnel, 917005 epoxy body ATC probe, 475 ml each of pH 4.01 buffer, pH 7.00 buffer, pH 10.01 buffer, and pH Electrode Storage Solution, 080113 DO probe service kit, electrode stand, 080017 DO cal sleeve.

DP-1116A 1787.40/ea





ORION® 4-STAR BENCHTOP pH/CONDUCTIVITY METER



- ★ Measures conductivity, resistivity, salinity, TDS, pH, and temperature
- ★ Bright LCD display with ready indication
- ★ Automatic or manual temperature compensation
- ★ Accepts either 4-electrode or 2-electrode cells (conductivity)
- ★ pH calibration: accepts from 1 to 5 points with either automatic or manual temp compensation
- ★ IP54 dustproof and splashproof
- ★ Advanced self-tests and self-diagnostics
- ★ 200-point datalogger downloadable to printer or computer
- ★ Three-year replacement warranty

pH

Range	-2.000 to 19.999
Resolution	0.1/0.01/0.001
Relative Accuracy	±0.002

Conductivity

Range	0.000 to 3000 mS/cm, auto resolution, with cell constant dependence
Resistivity	0.0001 to 100 Megohms
Salinity	0.01 to 80.0 ppt NaCl equivalent
TDS	0.0 to 19900
Reference Temperature	15°, 20° or 25°C (default), autoranging
Cell Constant	0.001 to 199.9
Temperature	-5 to 105°C

ORION® 4-STAR pH/CONDUCTIVITY BENCHTOP METER ONLY

MO-1117 1181.70/ea

CONDUCTIVITY/pH KIT

Includes meter, 9157BNMD pH/ATC electrode, 013005MD conductivity cell, one pack 5x60 ml 1413 µS/cm conductivity standard, 475 ml each of pH 4.01, 7.00 and 10.01 buffers, pH Electrode Storage Solution and Electrode stand.

MO-1117A 1820.70/ea



Conductivity + pH

Portable

YSI MODEL 63 pH/CONDUCTIVITY METER

Measure pH, conductivity, salinity and temperature—all with one probe! The YSI 63 incorporates all the pH features of the YSI 60 and the proven performance of other YSI conductivity meters.

As with the YSI 60, the meter is supplied with an attached cable of 10, 25, 50 or 100 feet. pH readings are automatically temperature compensated. pH resolution: 0.01; accuracy: ±0.1. Temperature resolution: 0.1°C; accuracy: ±0.1°C. Conductivity resolution: 0.1 µS to 0.1 mS, depending on range; accuracy: ±0.5%. Meter warranty: 2 years. Probe warranty: 1 year.

Cable Length	NCL #	Price/ea
10 feet	YSI-6310	1035.00
25 feet	YSI-6325	1095.00
50 feet	YSI-6350	1165.00



ORION® 4-STAR PORTABLE pH/CONDUCTIVITY METER

- ★ Measures conductivity, resistivity, salinity, TDS, pH, and temperature
- ★ Bright LCD display with ready indication
- ★ Automatic or manual temperature compensation
- ★ Accepts either 4-electrode or 2-electrode cells (conductivity)
- ★ pH calibration: accepts from 1 to 5 points with either automatic or manual temp compensation
- ★ IP54 dustproof and splashproof
- ★ Advanced self-tests and self-diagnostics
- ★ 200-point datalogger downloadable to printer or computer
- ★ Three-year replacement warranty



Conductivity	
Range	0.000 to 3000 mS/cm
Resistivity	0.0001 to 100 Megohms
Salinity	0.01 to 80.0 ppt NaCl equiv.
TDS	0 to 19999
Reference Temp.	15°, 20° or 25°C (default)
Cell Constant	0.001 to 199.9
Temperature	-5 to 105°C

pH	
Range	-2.000 to 19.999 or 0 to 14
Resolution	0.1/0.01/0.001
Relative Accuracy	±0.002

ORION® PORTABLE pH/CONDUCTIVITY METER ONLY

MO-1217 823.50/ea

pH/CONDUCTIVITY KIT

Includes meter, 9107BNMD pH/ATC electrode, 013005MD 4-electrode conductivity cell, one 10-pack each of pH 4.01, pH 7.00, and pH 10.01 buffers, one pack 5 x 60 ml 1413 µS/cm conductivity standard, 5 x 60 ml pH electrode storage solution, one 50 ml beaker, and hard field case.

MO-1217A 1250.10/ea

ORION® 3-STAR PORTABLE pH METER

- ★ Measures pH, mv, and Temp
- ★ Bright LCD display with ready indication
- ★ Automatic or manual calibration (1 to 5 points)
- ★ Automatic buffer recognition
- ★ IP67 waterproof (may be immersed for up to one hour with no water incursion)
- ★ Advanced self-tests and self-diagnostics on meter
- ★ Convenient Quick-Start Guide is located on the back of the meter
- ★ 200 point datalogger which is downloadable to printer or computer via an optional RS232 or USB interface accessory
- ★ Three-year replacement warranty
- ★ One input (single channel)
- ★ BNC pH electrode connector
- ★ Mini-DIN connector for temp probe

pH Range	-2.000 to 19.999
Resolution	0.1 / 0.01 / 0.001
Relative Accuracy	±0.002
Temperature	-5 to 105°C
mV Range/Resolution	±1999.9 / 0.1

ORION® 3-STAR PORTABLE pH METER ONLY**PO-1212 462.65/ea****PORTABLE METER pH APPLICATION KIT**

Includes the meter, 9107BNMD triode, one 10-pack each of pH 4.01 buffer, pH 7.00 buffer, pH 10.01 buffer, and pH Electrode Storage cell, one 50 ml beaker, and hard field case.

PO-1212A 661.50/ea

PORTABLE WATERPROOF pH METER WITH 200-POINT DATALOGGER

Portable waterproof pH/mv/Temp meter developed for hard on-site work. Features a **200-point datalogger**, large LCD display, RS232 interface, automatic temperature compensation, and autorecognition of 4, 7, and 10 buffers, using a 2-point calibration. Battery operated using four "AA" disposable batteries. Datalogging intervals are between 5 seconds and one hour. Logged data can be viewed from meter display or can be downloaded to a PC using optional software and cable. Meter has DIN connector which is compatible with WTW probes. Meter can accommodate any BNC probes using a BNC to DIN adapter (not included). Probe is gel-filled and submersible.

3 Year warranty on meter.



Description	NCL #	Price/ea
Meter, Sentex 41 pH/Temp probe with 3 ft cable, case, electrode support, buffers (4, 7, 10), storage solution, and beaker	WT-794	834.00
Replacement Sentex 41 pH/Temp probe with 3 ft cable	WT-635	233.00
Replacement Sentex 41 pH/Temp probe with 10 ft cable	WT-636	272.00
BNC to DIN adapter	WT-BUC1	76.25
WTW data collection software and cable	WT-852	175.00

DIGITAL pH/TEMP METER FOR PLATING BATHS

Digital waterproof meter has stainless steel probe (not submersible) with 3 ft cable. Shows pH and Temp simultaneously. pH range -2.00 to 16.00, Temp range -5.0 to 105.0°C (23 to 221°F).

Meter can be switched between °C and °F. pH resolution is 0.01 and accuracy is ±0.02. Temp resolution is 0.1° and accuracy is ±0.5° up to 60°C and ±1° above that. Automatic temperature compensation. Features automatic calibration using 2 buffers (any two of 4.00, 7.00, and 10.00). Powered by three 1.5V batteries. Package includes meter, probe, batteries, instruction manual, and carrying case. DIN Connector.

PH-131 **373.50/ea**

REPLACEMENT PROBE FOR ABOVE

PH-911 **120.60/ea**

SHOCK-PROOF BOOT

PH-007 **31.50/ea**





YSI MODEL 60 pH METER

The new YSI 60 pH system carries a waterproof rating of IP65, and its pH electrode and cable are completely submersible. The combination pH electrode uses a waterproof BNC connector, making it easy to replace in the field. The YSI 60 stores up to 50 sets of data and features a large backlit LCD, which is easy to read in any light. pH readings are automatically temperature compensated.

pH resolution: 0.01; accuracy: ± 0.1

Temperature resolution: 0.1°C; accuracy $\pm 0.1^\circ\text{C}$

Meter warranty: 2 years

Probe warranty: 1 year

Cable Length	NCL #	Price/ea
10 feet	YSI-6010	706.50
25 feet	YSI-6025	760.50
50 feet	YSI-6050	787.50

REPLACEMENT pH ELECTRODE
YSI-131133 195.00/ea

ECONOMY PORTABLE METERS FROM YSI



YSI pH/mv/TEMP METER, ECONOMY MODEL

The pH100 is YSI's new economy version of a portable pH/mv/Temp meter. It features a watertight housing, internal buffer recognition, auto/manual temperature compensation, low battery indicator, and auto hold of readings. Resolution is 0.01 in the pH mode and 1 mv in the millivolt mode. The instrument has a one year warranty. Probes are available for purchase separately.

YSI-100 (meter only) 250.00/ea

PROBES FOR THE YSI-100

Description	NCL #	Price/ea
Waterproof pH/Temp probe with 3 ft cable	DY-1001	94.05
Waterproof pH/Temp probe with 12 ft cable	DY-1004	118.75



SOFT-SIDED CARRYING CASE

There is a velcro-grip strap to hold the meter in place, an open area for cables, membrane kits, and elastic compartments for miscellaneous.

YSI-485 26.00/ea

POCKET pH TESTERS

BARGAIN HUNTERS NOTE: As with most other things in life, you get what you pay for when it comes to lab equipment. Pocket pH testers do NOT take the place of a pH meter. They are useful for screening of samples only. They SHOULD NOT be relied upon for "reportable results". Here we have listed two of the more reliable ones from a standpoint of accuracy and reliability.

**pHep®4 & pHep®5**

Features:

- ✓ Waterproof, can be fully immersed in water
- ✓ Display pH/Temp at the same time
- ✓ Replaceable pH electrode
- ✓ Stability indicator
- ✓ Auto shut-off
- ✓ Automatic 1 or 2 point calibration using 4, 7, and 10 buffers
- ✓ Automatic temperature compensation
- ✓ Supplied with protective cap, electrode removal tool, batteries, and instructions

	pHep®4	pHep®5
NCL #	PH-127	PH-128
pH Range	0.0 to 14.0	0.00 to 14.00
pH Resolution	0.1	0.01
pH Accuracy	±0.1	±0.05
Price	76.50/ea	84.60/ea

REPLACEMENT pH SENSOR FOR THE ABOVE
PH-128A 41.40/ea

PORTABLE pH/ORP/ISE METER

- ◆ Rugged and waterproof
- ◆ pH: Accepts up to 3 calibration points with autorecognition of 5 standard buffers
- ◆ Displays pH and Temp at the same time
- ◆ ISE: Accepts up to 5 calibration points.

Parameter	Range	Resolution	Accuracy
pH	-2.00 to 16.00	0.1/0/01	±0.01
mv	-2000 to 2000	0.1	±0.2
Temperature	-20.0 to 120.0°C	0.1°C	±0.4°C
ISE	0.001 to 1000 ppm	0.01, 0.1, 1	±0.5% for monovalent ions ±1% for divalent ions



Works with: ammonia, bromide, calcium, carbon dioxide, cupric fluoride, iodide, lead, nitrate, potassium, silver, sulfate, and sulfide probes.

Meter package includes meter, pH probe, ATC probe, rechargeable batteries, inductive charger, and rugged carrying case.

Meter package **HI-98184** **855.00/ea**
Replacement pH/ATC probe **HI-72911** **131.40/ea**

**ECONOMY HAND-HELD pH/mv/TEMP METER**

- ◆ Heavy duty waterproof design
- ◆ Large LCD readout in a low-profile case
- ◆ Displays pH and Temp at the same time
- ◆ Automatic 1 or 2 point calibration out of 5 stored buffers (4.01, 6.86, 7.01, 9.18, and 10.01)
- ◆ ±0.01 resolution

Meter package includes meter, pH probe, ATC probe, pH 4 & 7 buffer solutions, 100 ml plastic beaker, manual and rugged carrying case.

Meter package **PH-9125** **337.50/ea**
Probe **HI-1230B** **76.50/ea**
Temp Probe **HI-76692W** **56.70/ea**





ORION® 2-STAR BENCHTOP pH METER

- ★ Measures pH, mv, and Temp
- ★ Bright LCD display with ready indication
- ★ Automatic or manual calibration (1 to 3 points)
- ★ Automatic buffer recognition
- ★ IP54 dustproof and splashproof
- ★ Three-year replacement warranty


ORION® 2-STAR pH BENCHTOP METER ONLY

PO-1111 445.55/ea

APPLICATION KIT:

Includes the meter and 9107BNMD epoxy gel-filled triode

PO-1111G 528.30/ea

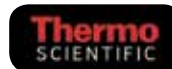

ORION® 3-STAR BENCHTOP pH METER

- ★ Measures pH, mv, and Temp
- ★ Bright LCD display with ready indication
- ★ Automatic or manual calibration (1 to 5 points)
- ★ Automatic buffer recognition
- ★ IP54 dustproof and splashproof
- ★ Advanced self-tests and self-diagnostics on meter
- ★ 200 point datalogger which is downloadable to printer or computer via an optional RS232 or USB interface accessory
- ★ Three-year replacement warranty

ORION® 3-STAR pH BENCHTOP METER ONLY

PO-1112 552.60/ea

APPLICATION KIT:

Includes the meter, 9157BNMD triode, electrode stand, 475 ml each of pH 4.01 buffer, pH 7.00 buffer, pH 10.01 buffer, and pH electrode storage solution.

PO-1112A 799.20/ea

ECONOMY BENCHTOP pH/mv/TEMP METER

- ◆ Large, easy to read LCD display
- ◆ Displays pH and Temp at the same time
- ◆ Automatic 1 or 2 point calibration
- ◆ ±0.01 resolution
- ◆ Meter package includes refillable pH probe, one 30 ml bottle of filling solution, temperature probe, and probe holder

Description	NCL #	Price/ea
Meter Package	HI-2211	351.00
pH Probe	HI-1131B	99.90
Temp Probe	HI-76692W	56.70
Probe Holder	HI-76404	63.90



LABORATORY pH/ISE METERS WITH TWO INPUTS**ORION® DUAL STAR METER**

2-channel (2 input) pH/ISE meter

Use one channel for pH and the other for ISE measurements, or both channels for ISE measurements.



- ★ Displays the results on both channels simultaneously
- ★ USB and RS232 data ports
- ★ 1000 point datalogger with time & date stamp
- ★ Calibration graph with slope data for viewing or printing
- ★ Calibrations with up to 5 points possible
- ★ Three-year warranty

**DUAL STAR pH/ISE METER ONLY**

PO-2115 1255.50/ea

DUAL STAR METER WITH ROSS ULTRA

Includes meter, 8102BNWP probe, electrode stand, 927007MD ATC probe, and RS232 cable.

PO-2115A 1466.10/ea

	pH	mV	ISE
Range	-2 to 19.999	-1999 to 1999	0 to 19999
Resolution	0.1, 0.01, 0.001	0.1	1, 2, or 3 significant digits
Relative Accuracy	±0.002	±0.2	0.05%

DUAL STAR pH/ISE METER KIT

Includes meter, 8102BNWP probe, electrode stand, 927007MD ATC probe, and RS232 cable. Also includes one pint each of pH 4, pH 7, pH 10 buffers, and pH electrode storage solution.

PO-2115B 1509.30/ea



DUAL CHANNEL pH/ISE METER

Research grade laboratory benchtop meter features a 240 x 320 dot-matrix color display. The large LCD screen allows easy navigation through menus and options without having to consult the manual. For pH, the meter features a 5 point calibration with a choice of custom or standard buffers and provides the user with the exclusive Calibration Check™ diagnostic system. For ISE measurements, the meter allows direct calibration and measurement. Meter accepts up to 5 point calibration; choose from 5 fixed standard solutions and 5 user defined standards.

- ◆ Displays one or both channels
- ◆ Automatic endpoint detection
- ◆ Displays offset and slope
- ◆ 50,000 point datalogger
- ◆ Datalog on command, or by user-set intervals
- ◆ PC Connection Opto-isolated USB and RS232
- ◆ Manual or automatic temp compensation



Meter package includes meter, glass body refillable pH electrode, temperature probe, AC adapter, electrode holder, and probe filling solution.

HI-4222

1170.00/ea

	pH	mv	Temperature	ISE
Range	-2.00 to 16.00	-2000 to 2000	-20.0 to 120.0°C	1×10^{-7} to 9.99×10^{10}
Resolution	0.1 / 0.01 / 0.001	0.1	0.1°C	1, 0.1, 0.01
Rel. accuracy	±0.1 / 0.01 / 0.002	±0.2	±0.2°C	±0.5%, ±1.0%



ORION® 4-STAR PORTABLE pH/ISE METER

- ★ Measures pH, mv, and Temp
- ★ Bright LCD display with ready indication
- ★ Accepts ISFET electrodes
- ★ Automatic or manual calibration (1 to 5 points)
- ★ Automatic buffer recognition
- ★ Reads concentration directly
- ★ Stores up to 10 methods with password protection
- ★ Point-to-point ISE calibration
- ★ Automatic blank correction
- ★ Advanced self-tests and self-diagnostics on meter
- ★ 200 point datalogger which is downloadable to printer or computer via an optional RS232 or USB interface accessory
- ★ IP67 waterproof
- ★ Three-year replacement warranty

PORTABLE METER pH APPLICATION KIT

Includes the meter, 9107BNMD triode, one 10-pack each of pH 4.01 buffer, pH 7.00 buffer, pH 10.01 buffer, one pack of 5 x 60 ml pH electrode storage solution, and hard field case.

PO-1215A 1001.70/ea

pH

Specifications

Range	-2.000 to 19.999 or 0 to 14 in aqueous solutions
Resolution	0.1 / 0.01 / 0.001
Relative Accuracy	±0.002

ORION® 4-STAR PORTABLE pH/ISE METER ONLY

PO-1215 830.70/ea

ISE

Range	0 to 19900
Resolution	1, 2, or 3 significant digits
Relative Accuracy	±0.2mV or 0.05%, whichever is greater
Units	M, mg/l, %, ppb, or no units
mv/Rmv	±1999.9
Temperature	-5 to 105°C





ORION® 4-STAR BENCHTOP pH/ISE METER

- ★ Measures pH, mv, and Temp
- ★ Bright LCD display with ready indication
- ★ Accepts ISFET electrodes
- ★ Automatic or manual calibration (1 to 5 points)
- ★ Automatic buffer recognition
- ★ Reads concentration directly
- ★ Stores up to 10 methods with password protection
- ★ Point-to-point ISE calibration
- ★ Automatic blank correction
- ★ Advanced self-tests and self-diagnostics on meter
- ★ 200 point datalogger which is downloadable to printer or computer via an optional RS232 or USB interface accessory
- ★ IP54 dustproof and splashproof
- ★ Three-year replacement warranty

**ORION® 4-STAR pH/ISE BENCHTOP
METER ONLY**
PO-1115 1075.50/ea

Specifications

pH	Range	-2.000 to 19.999 or 0 to 14 in aqueous solutions
	Resolution	0.1 / 0.01 / 0.001
	Relative accuracy	± 0.002
ISE	Range	0 to 19900
	Resolution	1, 2, or 3 significant digits
	Relative accuracy	± 0.2mV or 0.05%, whichever is greater
mV/ORP	Units	M, mg/l, %, ppb, or no units
	Resolution	0.1
	Relative accuracy	± 0.2mV or 0.05%, whichever is greater
	mv/Rmv/EH ORP	± 1999.9
	Temperature	-5 to 105°C

ENVIRONMENTAL APPLICATIONS KITS**pH**

Kit includes: meter, 8102BNUWP ROSS Ultra™ pH electrode, 917005 epoxy ATC probe, electrode stand, 475 ml each of pH 4.01 buffer, pH 7.00 buffer, pH 10.01 buffer, and pH electrode storage solution.

KO-5001 1328.40/ea

AMMONIA IN WASTE WATER

Kit includes meter, 9512BNWP ammonia electrode (no reference probe required), one 60 ml bottle of O-951202 probe filling solution, electrode stand, 096019 stirrer and paddle, one pack O-951204 (20 loose membranes), 475 ml each Ammonia ISA Buffer, and 100 ppm Standard.

KO-109 1971.00/ea

FLUORIDE IN DRINKING WATER

Kit includes meter, 9609BNWP ion-plus fluoride electrode (no reference probe required), one 60 ml bottle of O-900061 probe filling solution, electrode stand, 096019 stirrer and paddle, 475 ml each of 1 ppm Fluoride Standard w/TISAB II, 2 ppm with TISAB II, 10 ppm with TISAB II, and one gallon of TISAB II.

KO-103 2132.10/ea

