

## EXCELLENT QUALITY "ECONOMY" FECAL WATER BATHS




<b>NCL #</b>	<b>BP-101A</b>
<b>Reservoir Capacity</b>	28 liters
<b>Inside Dimensions</b>	10.25" L x 11.75" W x 8" D
<b>Outside Dimensions</b>	14" L x 22" W x 13.5" H
<b>Temp Range</b>	Amb+5 to 150, preset@44.5°C
<b>Temp Uniformity</b>	±0.05°C
<b>Controller</b>	microprocessor
<b>Electrical</b>	120V, 8.3 amps
<b>Shipping Weight</b>	49 lbs
<b>Price</b>	<b>1795.00/ea</b>

## PRECISION WATER BATHS

These baths meet APHA, AWWA, EPA, and WEF specifications, and come factory calibrated to 44.5°C, but can easily be re-adjusted by the user if necessary. All stainless steel construction with epoxy-powder paint coated exterior surfaces. Stainless steel gable covers included with both models. The smaller bath has an external 50 gph pump. The larger bath has a pump inside the housing. Both baths have a temperature sensitivity of ± 0.05°C and a temperature uniformity of ±0.2°C when used with the gable cover provided.

**Racks are not included.**

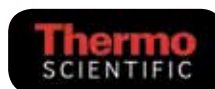

<b>NCL #</b>	<b>BP-100</b>	<b>BP-150</b>
<b>Capacity (gallons)</b>	4.8	9.1
<b>Inside Dimensions</b>	12" W x 14" D x 8" H	27" W x 12" D x 7 1/2" H
<b>Outside Dimensions</b>	15 1/2" W x 18 1/2" D x 9 1/2" H	36" W x 14" D x 9 1/2" H
<b>Power Requirements</b>	110VAC, 665 W	110VAC, 1550 W
<b>Net Weight (empty)</b>	27 lbs	46 lbs
<b>Price</b>	<b>2268.00/ea</b>	<b>2520.00/ea</b>

## REPLACEMENT PUMP FOR BP-100 WATER BATH BP-101 799.00/ea



# WATER BATH ACCESSORIES

## COLIFORM THERMOMETER FOR WATER BATHS

Range is 34.5 to 46°C in 0.1°C divisions.  
Total length is 237 mm with 76 mm immersion. Critical temperatures of 35, 37, 44.5, and 45.5°C are marked with arrows. Blue liquid-filled.



**TC-35A**                      **68.60/ea**

## CERTIFIED COLIFORM THERMOMETER FOR WATER BATHS

Same as the above, except certified @ 44.5°C.



**TC-35CERT**                      **103.60/ea**

## STAINLESS STEEL PETRI DISH RACK FOR WATER BATHS

Holds 96-50 mm petri dishes.

**BP-103**                      **222.00/ea**



## WATER BATH CONDITIONER

Aqua-Clear Water Conditioner stabilizes water in water baths and other water systems. Just 2 ml per liter of water provides long lasting protection against the growth of algae, bacteria, and molds which often contaminate standing water. It also protects pumps and thermostats and helps prevent tubing from clogging. Aqua-Clear Water Conditioner adds a blue tint to the water which will fade after about 3 months reminding you to change the water and add new conditioner. 100 ml polyethylene bottle with cap.

**DB-179**                      **66.00/ea**



# COLIFORM DRY BATHS FOR FECAL AND TOTAL COLIFORMS

## MILLIPORE DRY BATH FOR COLIFORM DETERMINATIONS

Holds up to 48 petri dishes and does not require special petri dishes. Designed for laboratory incubation of 47 mm microbiological test filters at 35 or  $44.5 \pm 0.2^\circ\text{C}$ . Temperature is thermostatically controlled and monitored by a thermometer which is included. Choice of five set temperatures: 30, 35, 37, 41, and  $44.5^\circ\text{C}$ . Molded polyethylene case, coated aluminum chamber and shelves, foam insulation, and clear plastic inner door. External dimensions: 13" H x 14.6" W x 11.7" D. 115VAC.

**IM-101**                      **1885.00/ea**

## REPLACEMENT THERMOMETER

Teflon-coated thermometer has a range of 28.0 to  $45.0^\circ\text{C}$  in  $0.2^\circ\text{C}$  increments. Equipped with two plastic clips specifically designed to attach to shelf of Millipore Dry Baths. Overall length: 160 mm.

**NOTE:** This thermometer contains mercury.

**IM-119T**                      **96.35/ea**

## REPLACEMENT DOORKNOB

**IM-101DK**                      **10.00/ea**

